
IMPIANTO DENTALE GOLD • GOLD DENTAL IMPLANT

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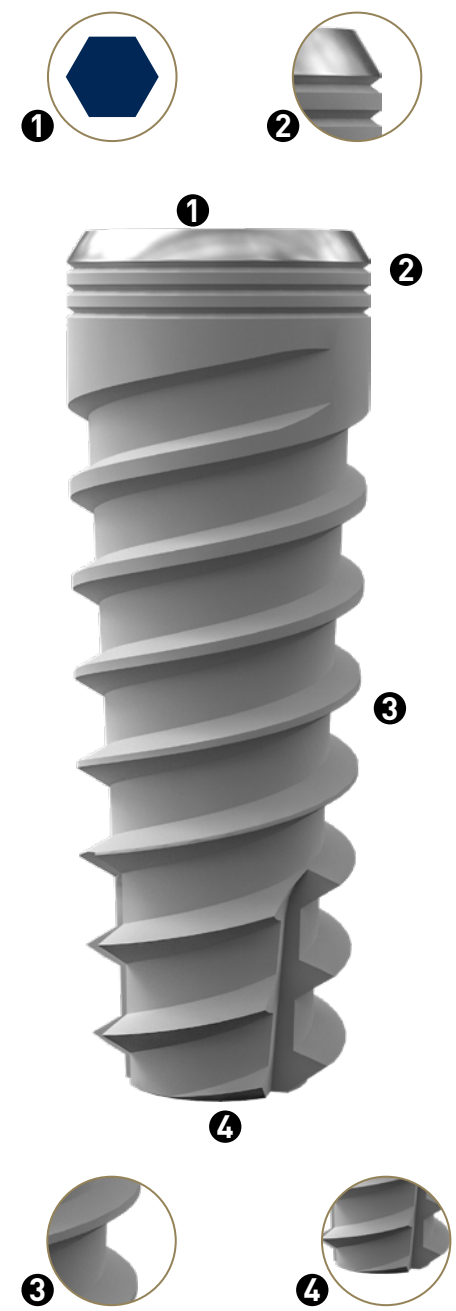GOLD is the implant line designed by Titanium Component to meet the needs of modern implantology, especially for immediate loading and post extraction, but also for anyone looking for a modern and versatile system that always ensures high primary stability, even in the most complex situations. The GOLD line is also characterized by a simple surgical protocol, and a modular surgical kit based on the specific clinical needs. Thanks to its technical characteristics, GOLD line implant provides many clinical advantages compared to conventional implants:

- The high conicity and the shape of the spiral ensure a high self-tapping effect, which results in a high bone compaction and an excellent primary stability. These advantages are particularly evident in soft bone, where GOLD line implants allow fast and secure prosthetic restorations.
- The implant does not require tapping, minimizing surgical stages and allowing, at the same time, the maximum preservation of the cortical bone, with minimal trauma even in complex rehabilitation cases.
- Thanks to the absence of machined collar, GOLD line implants can also be inserted over a millimeter below the bone crest, thereby improving the stability and limiting the occurrence of peri-implantitis.
- The spiral, more accentuated at the apex, ensures a very good grip and a high directionality even in post extraction, limiting the need for bone regeneration interventions.

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Gold line implants are packaged in a dual ampoule, inserted in a specific plastic blister for Gamma rays sterilization. In this way, Titanium Component provides maximum resistance of sterility, while ensuring patient health. The attention to sterility maintenance, however, is useless if not accompanied by a safe surgical procedure; in this regard, we have realized the cap of the inner ampoule in AISI630 medical steel, with a hexagonal shape designed to engage perfectly on the Titanium Component ratchet. In this way the assembly of the implant mounter on the ratchet can be carried out when the ampoule is still closed, completely eliminating the need to manipulate the implant by the operator.

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
TITANIUM COMPONENT sas
Via Socrate, 26 • 20128 • Milano • Italy
Telephone (+39) 02 49.75.06.34
mail: titaniumcomponent@gmail.com
web: www.protesicacompatibile.com

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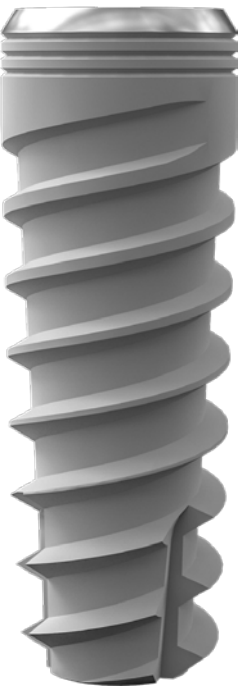KIT FRESE CON STOP • DRILL STOP KIT

The GOLD line has a wide range of prosthetic solutions with connection compatible with Zimmer® Screw Vent®. The connection is the same for all implant diameters, allowing to obtain, on larger diameters, a switching platform to further ensuring the effectiveness and durability of the prosthesis. An optimum result always comes from the attention to detail. For this reason, Titanium Component offers dental technicians a complete line of five different analogues, which reflect exactly the characteristics of the emergencies of each implant diameter, thus allowing maximum precision in the realization of individual crowns.

ANALOGO PER MUA IMPLANT REPLICA FOR MUA			LOCATOR® LOCATOR®				
Ø Unico/Unique Ø Unico/Unique	Angol./Angled Dritto/Straight	PL-AN0003 PL-AN0004		Ø Tutti/All Ø Tutti/All Ø Tutti/All Ø Tutti/All Ø Tutti/All Ø Tutti/All	H 1 H 2 H 3 H 4 H 5 H 6	PL-LC0001/1 PL-LC0001/2 PL-LC0001/3 PL-LC0001/4 PL-LC0001/5 PL-LC0001/6	
CILINDRO PER MUA CYLINDER FOR MUA							
Ø Unico/Unique Ø Unico/Unique	Angol./Angled Dritto/Straight	PL-CM0001 PL-CM0002					
CALCINABILE PER MUA CALCINABLE FOR MUA							
Ø Unico/Unique Ø Unico/Unique	Angol./Angled Dritto/Straight	PL-CN0005 PL-CN0006		Ø Tutti/All Ø Tutti/All	H 1,5 H 3,0	PL-LC0003 PL-LC0004	

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GOLD implant surgical kit is designed for maximum simplicity and intuitive use and includes, in addition to the dynamometric ratchet, eight drills equipped with notches for depth measurement. Titanium Component also provides a range of screw drivers and inserters that can be ordered separately, so as to ensure a high level of customization of the instruments.

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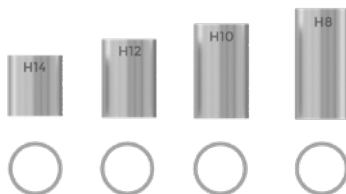
• It is recommended to insert the implant at least 0.5mm below the bone crest.
• In very soft bone it may be advisable to carry out the last step of drilling with the final drill taking care to maintain a depth of 2mm infer than the selected implant. In this manner, the bone is more compact in correspondence of the apex of the implant to the benefit of the primary stability.